




SYSTEMS BIOLOGY FOR DRUG DISCOVERY


- **1. Download <http://ftp.genego.com/dist/CytoscapeNetPlugin.jar>**
- **2. Place the CytoscapeNetPlugin.jar file in the 'plugins' folder within the Cytoscape installation directory. For example, C:\Program Files\Cytoscape_v2.5.1\plugins.**
- **3. Open Cytoscape**
- **4. The integration function is then accessible from the 'Plugins' menu in Cytoscape (see below).**
- **5. NOTE: When sending content to Cytoscape, the plugin will occupy a MetaCore license seat. The user will need to log out of MetaCore to release the seat for another user to use.**




- Select Import network from MetaCore from the Plugins menu

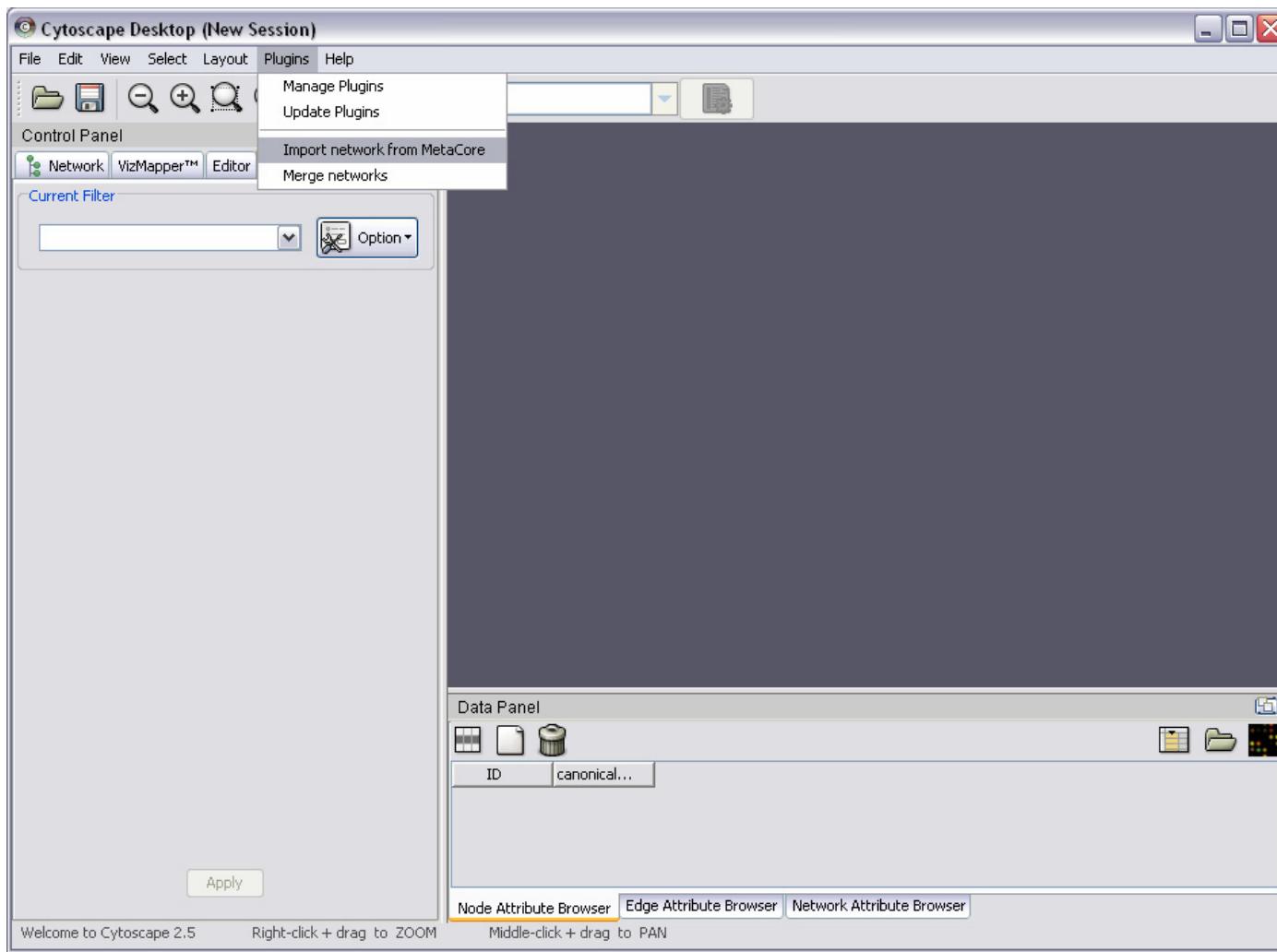

- You need to specify MetaCore server name, your login and password here

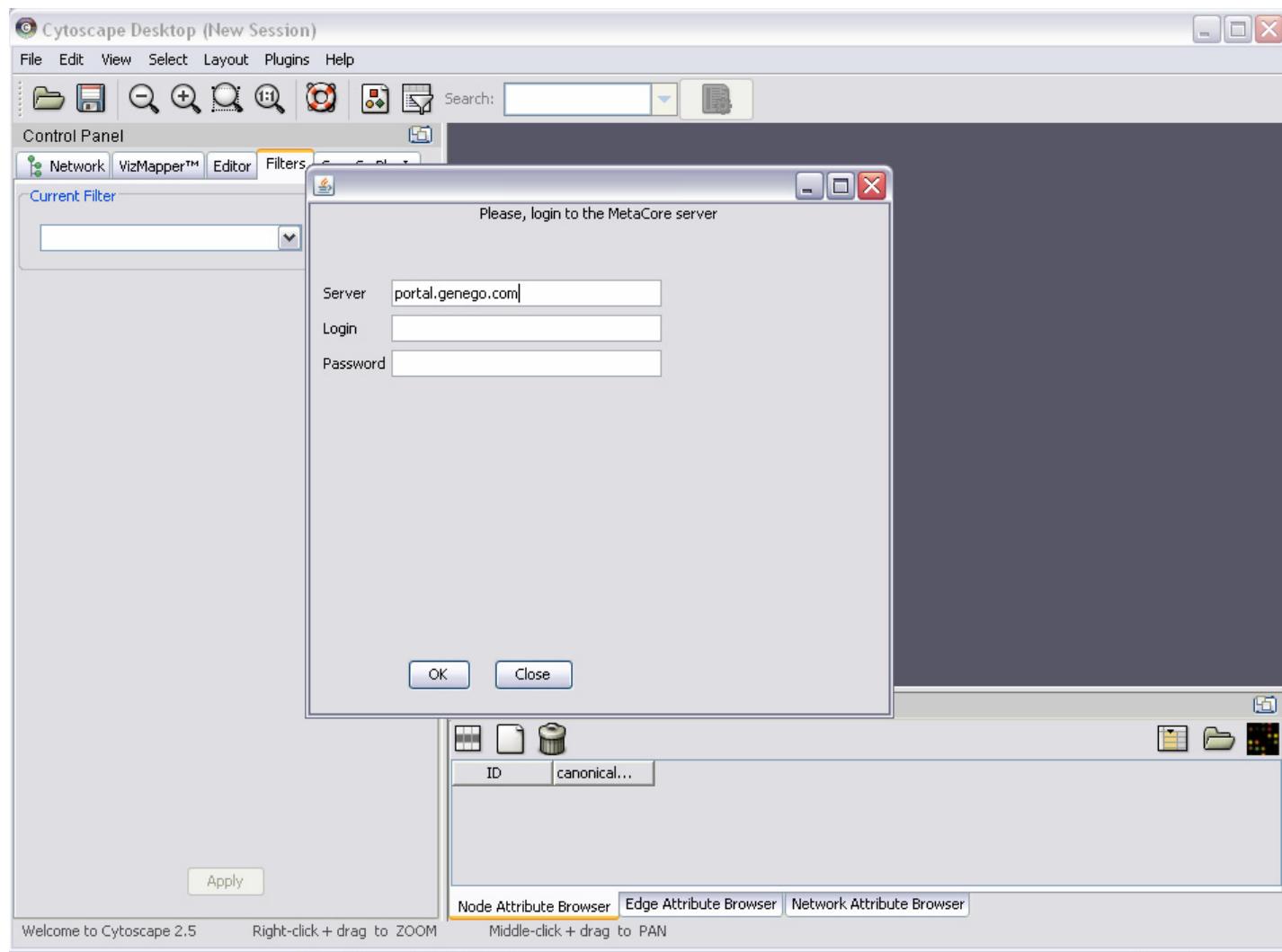

- Select network from your stored networks list

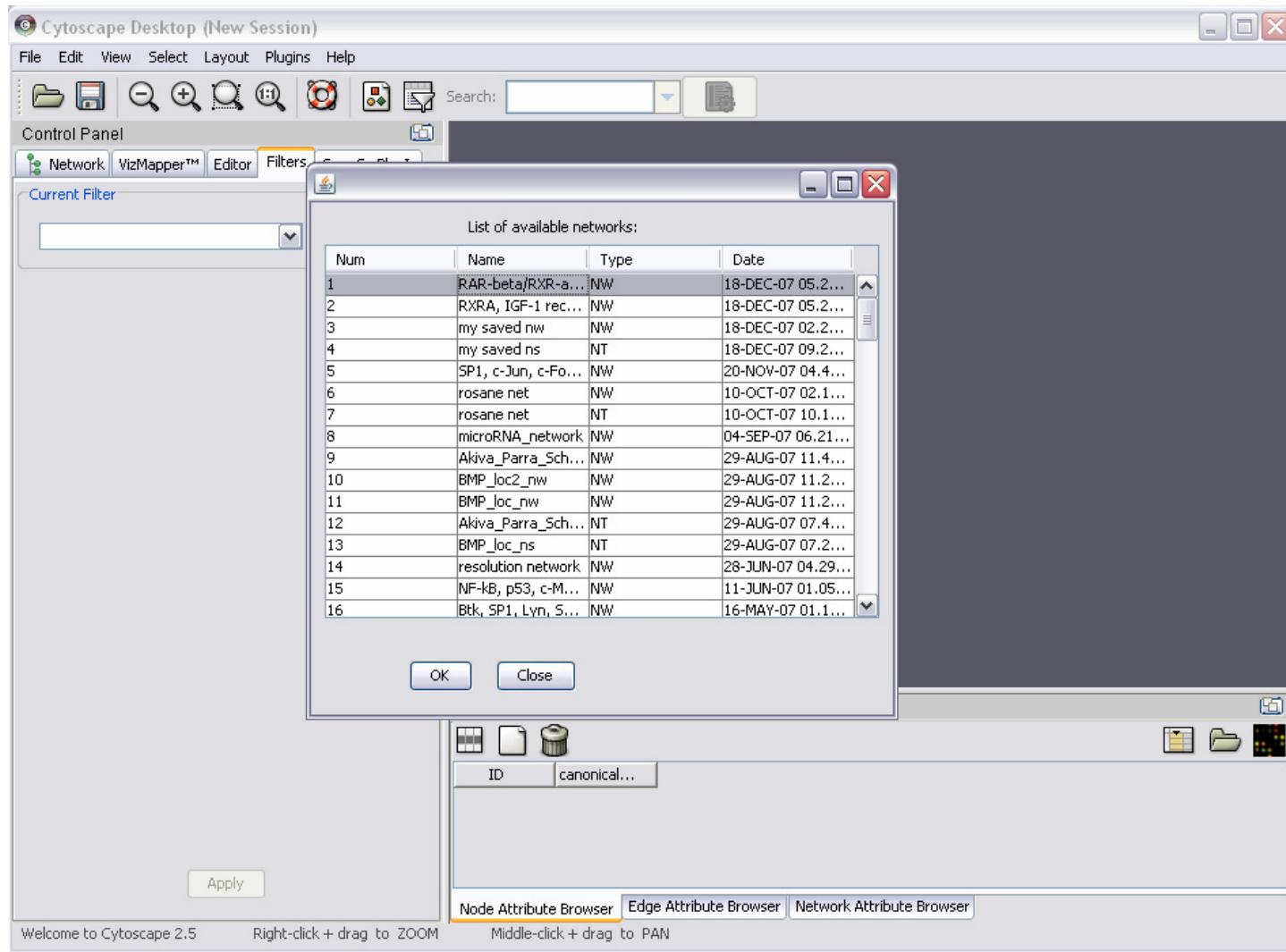


Welcome to Cytoscape 2.5


Right-click + drag to ZOOM

Middle-click + drag to PAN

- Selected network will appear in Cytoscape

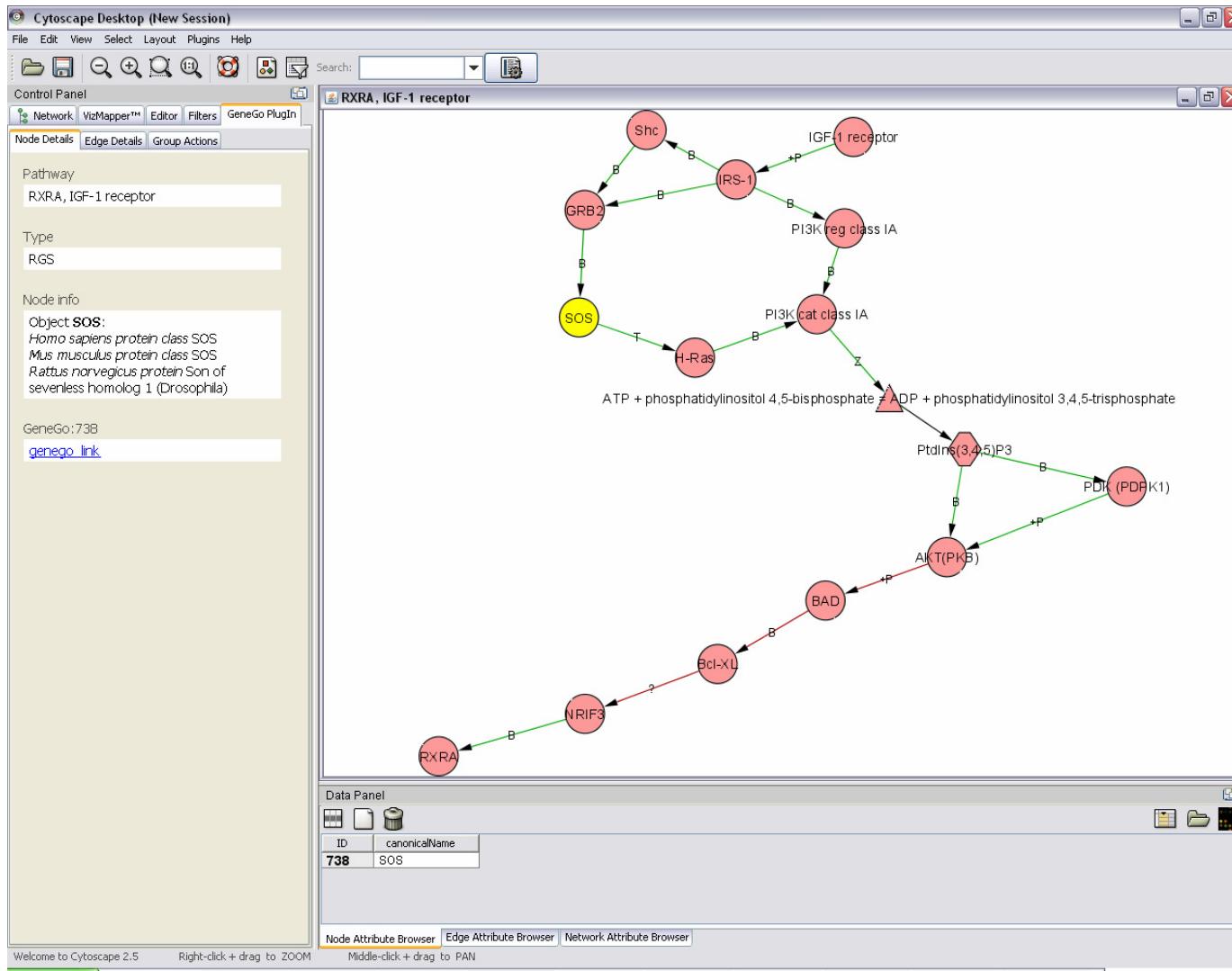

- More details for the selected node appear in the Node Details tab of the GeneGo Plugin tab. Click 'genego link' to open the MetaCore details page for the selected node

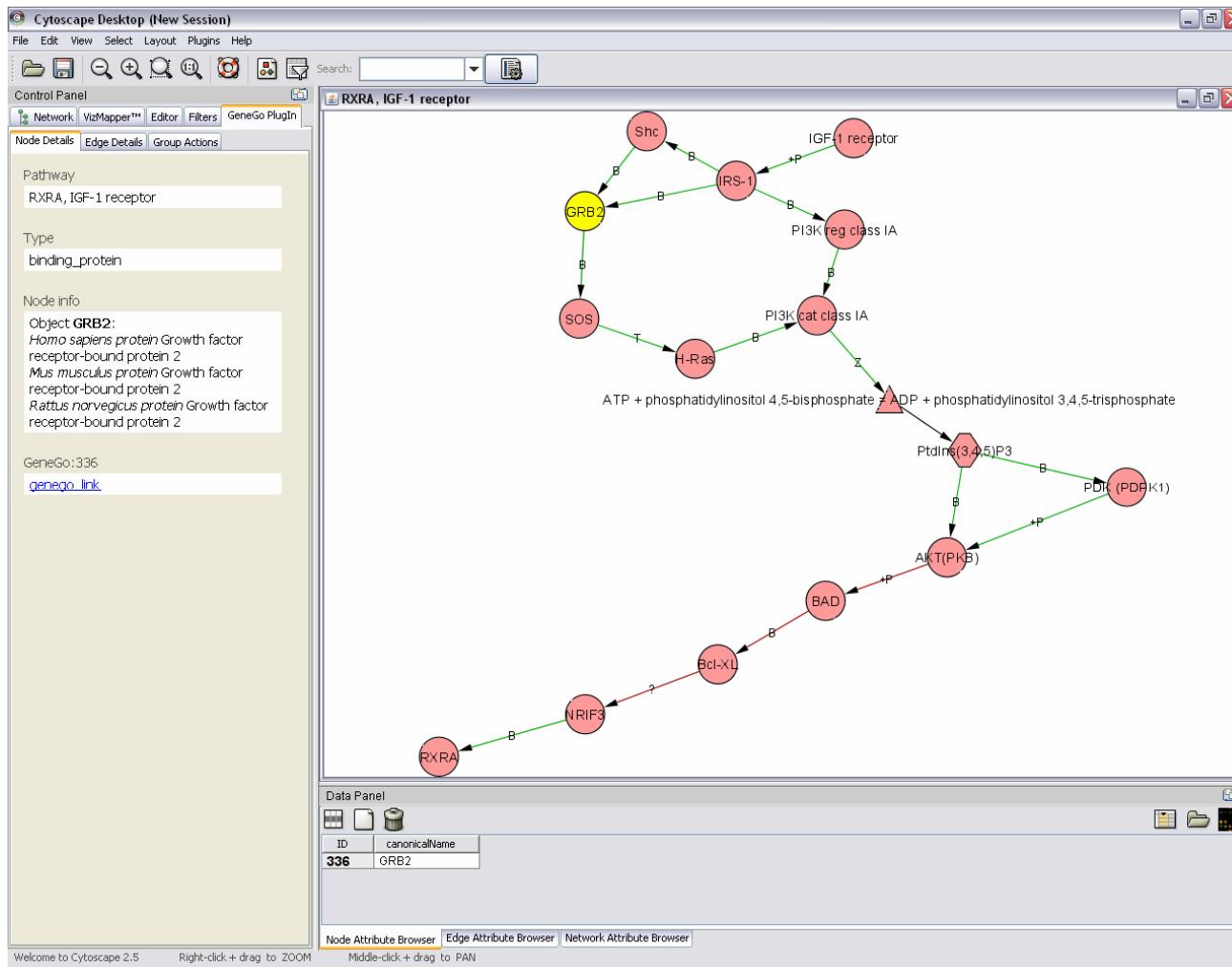

- More details for the selected edge appear in the Edge Details tab of the GeneGo Plugin tab. Click 'genego link' to open the MetaCore details page for the selected edge

